

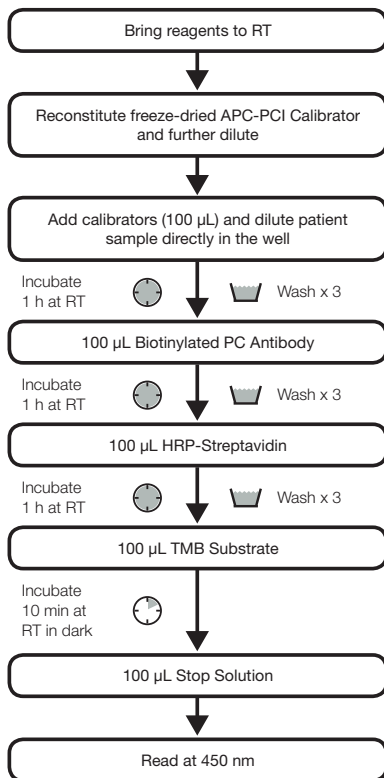


# Quick Guide

## APC-PCI ELISA Kit (KIT 040)

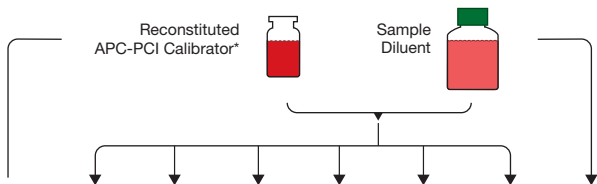

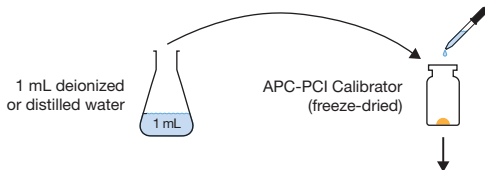

APC-PCI Sample Diluent	160 $\mu$ L 160 $\mu$ L	80 $\mu$ L 240 $\mu$ L	40 $\mu$ L 280 $\mu$ L	20 $\mu$ L 300 $\mu$ L	10 $\mu$ L 310 $\mu$ L	10 $\mu$ L 630 $\mu$ L	320 $\mu$ L
	Tube 1 1/2	Tube 2 1/4	Tube 3 1/8	Tube 4 1/16	Tube 5 1/32	Tube 6 1/64	Tube 7 -
	1 ng/mL	0.5 ng/mL	0.25 ng/mL	0.125 ng/mL	0.063 ng/mL	0.031 ng/mL	0 ng/mL
	100 $\mu$ L/ well x 2	100 $\mu$ L/ well x 2	100 $\mu$ L/ well x 2	100 $\mu$ L/ well x 2	100 $\mu$ L/ well x 2	100 $\mu$ L/ well x 2	100 $\mu$ L/ well x 2

Dilution Factor	1	2	4	8	16	32	64	-
APC-PCI	200 $\mu$ L	160 $\mu$ L	80 $\mu$ L	40 $\mu$ L	20 $\mu$ L	10 $\mu$ L	10 $\mu$ L	-
Sample Diluent	-	160 $\mu$ L	240 $\mu$ L	280 $\mu$ L	300 $\mu$ L	310 $\mu$ L	630 $\mu$ L	320 $\mu$ L

\*Reconstituted APC-PCI calibrator stock is stable for 24 h at 2-8°C.